


Identify an Unknown with Melting Point

Background

Melting point can be used to identify a substance and assess its purity. A pure substance has a well-defined melting point while an impure substance may have a lower, complex melting point that occurs over a broader range. The broader the melting range, the greater the impurity. Impurities disrupt crystal lattice structures and make it easier for a substance to transition from a solid to a liquid at a lower temperature. In this investigation, you will use a Wireless Melting Point Apparatus to identify a substance based on its melting point range. The Melting Point Apparatus heats up a sample in a slow and controlled manner to assure the temperature inside the sample tube is the same as the metal block heating it. The Heating Profile feature in SPARKvue is designed for accurate melting point range determination. The profile steps are visually represented in an image like the one below. The steps in the table are arranged in a typical sequence, however, steps can be reordered, added, or deleted as needed. The image is a model and the shape and completion time of the actual temperature plot will differ.

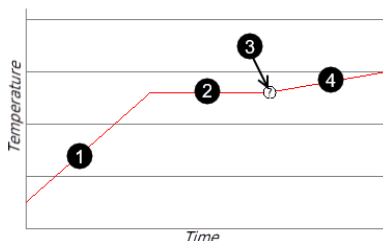

#	Step	Purpose
①	Quick Ramp	Rapidly approaches a specified temperature (in °C) in an asymptotic manner. Goal: Bring the sample's temperature a few degrees below the expected melting point.
②	Hold	Maintains the block at a steady temperature for a specified time in minutes. Goal: Allow the system and sample to reach thermal equilibrium. Choose a Quick Ramp temperature a few degrees below the melting point so the sample remains solid during Hold.
③	Prompt	Pauses the Heating Profile until the user clicks OK in the pop-up dialog that appears. Goal: Notifies the user when the sequence is about to go through the expected melting point so they can pay closer attention. If no action is taken within 30 minutes of the prompt's initial appearance, the sequence will automatically proceed to the next step.
④	Ramp	Slowly heat a sample to approach, meet, and exceed its anticipated melting point at a specified ending temperature and heating rate. Goal: Linearly approach and move through the expected melting point temperature at a rate between 1-2 °C per minute. Divide the increase in temperature over the time period to get the heating rate. Example: If Quick Ramp is set to (140 °C), and Ramp is set to (10 minutes, 150 °C), temperature will increase by 10 °C over 10 minutes therefore the heating rate equals 1 °C per minute.

Driving Questions

How can the melting point range of a substance help determine its identity? What are the signs that indicate when a phase change from solid to liquid first begins, and comes to completion?

Materials and Equipment

- Computer or mobile device
- Wireless Melting Point Apparatus
- USB Camera and camera adapter*
- Glass capillary tubes, one end closed, 1.5-mm OD
- Mortar and pestle
- Spatula
- Watch glass
- Sample(s) to be analyzed

*Optional; the camera adapter is included with the Wireless Melting Point Apparatus; check **pasco.com** to see which cameras are compatible with the apparatus, your system, and your software.

Safety

Follow these important safety precautions in addition to your regular classroom procedures:



- Wear gloves if available.
- Handle glass capillary tubes with care and dispose of used tubes as directed by your instructor.

Procedure

Part 1 - Setup

1. Begin your lab report with a 1-paragraph background description of the phenomena you will study in this investigation. Use your own words, and end with a clearly-stated question that can either be supported or denied with data collected during this investigation.
2. Plug the Wireless Melting Point Apparatus (MPA) into a wall outlet and switch the power on.
3. If available, connect a USB Camera to your computer or mobile device (before you open software).
4. Open SPARKvue and connect the MPA to your computer or mobile device, then open the Quick Start file called *Determining a Melting Point*.

NOTE: If the Quick Start file is not available, update your software to the latest version and start over.

5. For USB Camera users: Verify the camera is displaying a live preview on page 2 of the Quick Start file. Use the arrow at the top to navigate, click the Preview icon  in the right-side panel to turn on the camera feed, then return to page 1.
 - If the wrong camera is displayed, click the Preview icon to stop the live image feed. Open the Camera Tools menu . Go to the Settings gear and choose the correct camera. Click the Preview icon to show a live camera display, then return to Page 1.

Part 2 - Capillary Tube Preparation

1. Grind the sample to be analyzed with a mortar and pestle until the entire sample is a fine powder.
2. Transfer a small pile of the sample to the center of a watch glass, then tap the open end of the capillary tube into the sample as shown in Figure 1Ⓐ. Load it 2-3 mm deep.

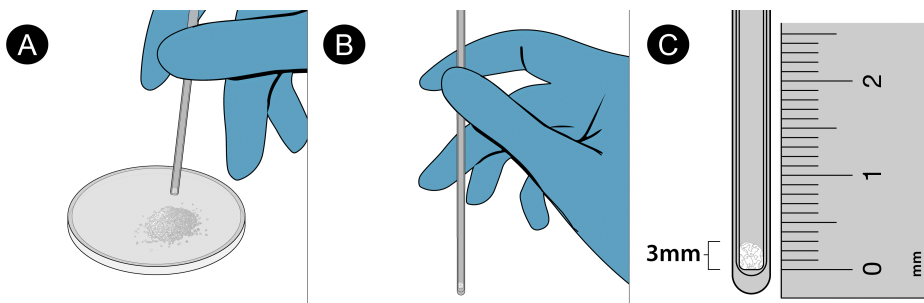


Figure 1. Collect a sample in a capillary tube

3. Invert the tube. Tap the bulbous tube end on the lab bench several times to compact the sample into the closed end (Figure 1Ⓑ).
4. Collect and compact additional sample if necessary to load the tube with a sample height of 2-3 mm (Figure 1Ⓒ). Do not exceed 3 mm.
5. Insert the closed end of the capillary tube into an open chamber in the MPA as shown in Figure 2. A magnified image of the sample will be visible on your screen if you are using a USB Camera (①) or in the eyepiece (②).

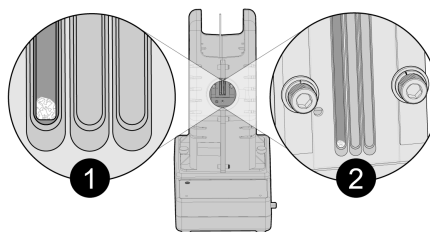





Figure 2. Melting Point Apparatus with sample in capillary tube

6. Choose either (a) or (b) depending on your viewing setup:
 - a. Eyepiece: Rotate the capillary tube and adjust the eyepiece position for the best focus and sample view. Remain on page 1 in the Quick Start file for the entire investigation.
 - b. USB Camera: If the video feed is not live, open the Camera Tools menu and turn on the Preview . Rotate the capillary tube and adjust the camera position to focus it for the best sample view, then return to page 1 in the experiment file.

Part 3 - Melting Profile and Melting Point Range

1. When working with an unknown substance, set up the Melting Profile to quickly estimate the melting point range:
 - a. Select the Ramp step near the bottom of the profile. Leave Time set as-is. Set Temperature to 350 °C.
 - b. Select Quick Ramp in the Melting Profile and delete it . Delete all remaining steps until only Ramp remains.
 - c. Locate the live temperature reading at the bottom-left corner. Assign lab partner roles: (1) sample monitor, and (2) temperature monitor. Lab partners should not switch roles during this investigation.
 - d. If using a camera, go to page 2 of the experiment file; otherwise remain on page 1.
 - e. Start collecting data. The sample monitor must let the temperature monitor know when melting begins and when it is complete. The temperature monitor will note the approximate temperatures at the start and end of melting.
 - f. Stop data collection a few degrees after melting is complete.
2. Select the Cool button above the Melting Profile. Allow the block to cool to at least 20 °C below the estimated first melting temperature, then turn off Cool.
3. Remove the used capillary tube from the MPA. Load a second capillary tube with the sample and set it in the same chamber. Adjust position and focus if needed.
4. Revise the profile on page 1 to apply a slower heating sequence based on the melting point range data determined in step 1:
 - a. Select the Select Profile Template button above the Melting Profile, then select Default from the drop-down list.
 - b. Select Quick Ramp in the Melting Profile. Change the temperature to a few degrees C below the observed temperature at which melting first began.
 - c. Set Hold to 3 minutes; allow Prompt to follow Hold.
 - d. Select Ramp. Set the ending temperature to a few degrees C above the observed temperature where melting ended, then determine the time required to set the rate of temperature increase to 1-2 °C per minute (see the Ramp step in the table in the Background section).
 - e. Use the arrows at the top of the display to navigate to experiment file pages. You can move back and forth among different pages at any time, even while collecting data. If using a camera, you may experience a brief delay in the image feed when you navigate between pages. Wait for the image to return. If the live image appears to freeze, open the Camera Tools menu  and turn on the Preview .
5. If using a camera, go to page 2; otherwise remain on page 1. Start collecting data; the active step in the Melting Profile will be highlighted. Monitor the sample and respond to the Prompt when it appears. If the sample has already begun to melt before the Prompt appears, stop collecting data. Cool the block as before and revise the Melting Profile as needed. Repeat steps 3 through 5 until your Heating Profile ensures the sample remains completely solid during Hold.

NOTE: Know the difference between glistening and melting. Melting begins when you can distinctly see some amount of liquid in the sample, especially on the sides of the capillary tube. Glistening involves light changing or sample settling, without the appearance of liquid. Melting is complete when the entire sample is in liquid form and is transparent.

6. You should now be in the Ramp step. Record the temperatures at which melting begins and ends, then allow the melting point profile sequence to complete (or, stop recording data when it is clear that melting is complete).
7. Repeat steps 2-6 until you have at least 2 runs with melting point range temperatures that agree within 0.50 °C.

Analysis

Answer all questions in your lab report.

1. Of the two runs that agree within 0.5 °C, select the lowest value for start of melting and the highest value for end of melting, then round to the nearest tenth degree C to determine the melting point range. What is the melting point range of your substance? Based on the melting point range, what is the likely identity of your substance?
2. Use the mean value in the melting point range to calculate your percent error. Show work.

Questions

Answer all questions in your lab report.

1. Comment on the purity of your sample based on your results, and discuss the validity of using percent error as calculated as a statement of sample purity.
2. Describe what you might see in the sample tube during heating if your unknown sample was composed of three different substances with distinctly different melting points.
3. Explain the importance of assigning roles in this investigation and maintaining those roles through the entire procedure.

Conclusion

Summarize your results to complete your lab report. Address the driving questions with generalizations from your data and observations. Use your understanding of intermolecular forces, colligative properties, phase changes, and data from this investigation to support your concluding statements. Comment on your accuracy based on your percent error calculations. Describe potential sources of error and their impacts on your results.